


# Neo Master Kit



This kit can be used to perform implant surgery for IS-III active, IS-II active, IT-III active.

## Compatible with IS-II active, IS-III active, IT-III active

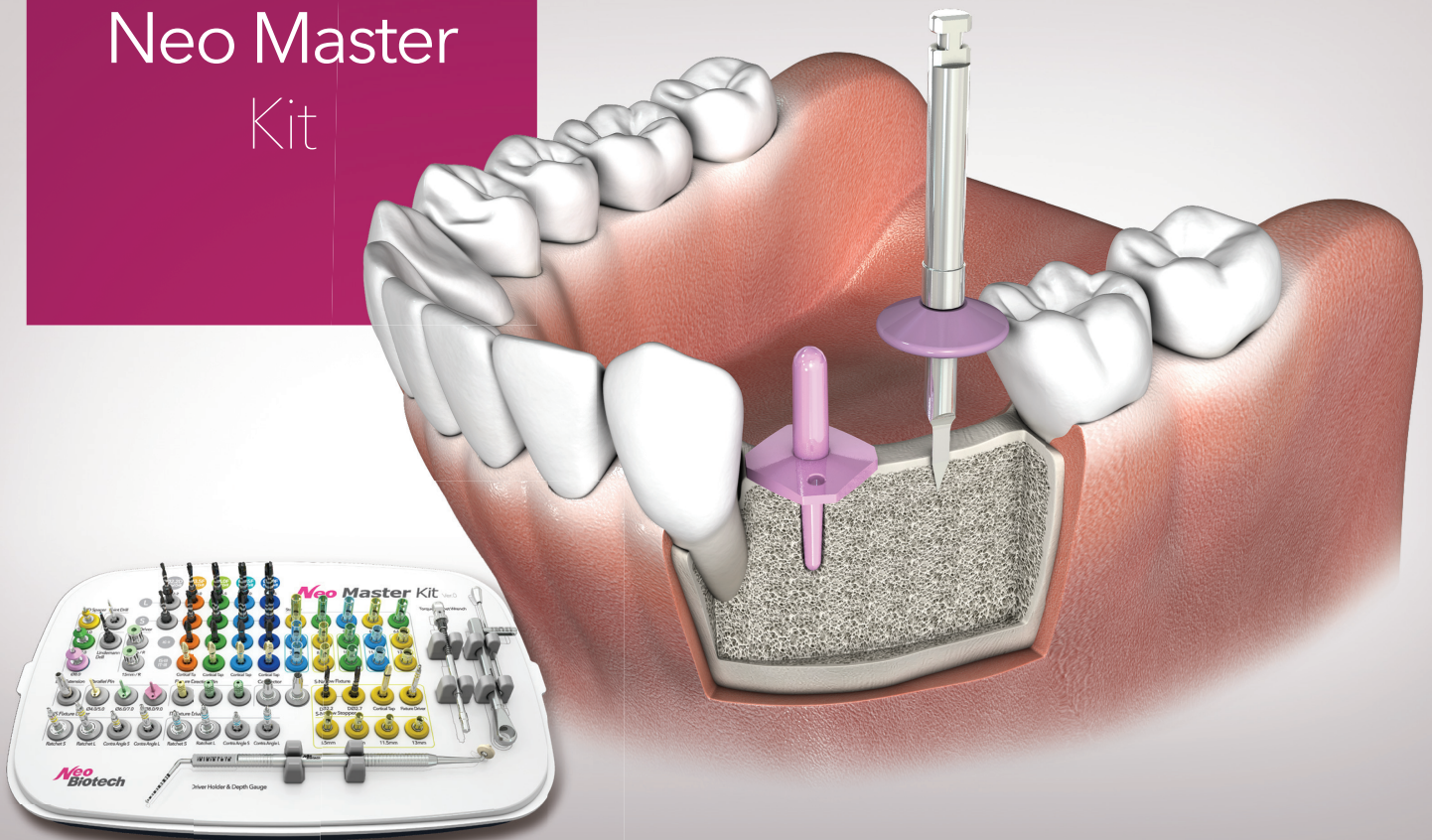


- Surgical drills and fixture drivers for fixture implantation are included.
- Applicable fixture types : IS-II active / IS-III active / IS-III active S-Narrow / IT-III active
  - For fixtures with  $\varnothing 5.5$  and more, S-Wide kit is required.

## Additional Drill Stopper




- Slots for frequently used drill stoppers are doubled;
- Additional drill stoppers are separately purchasable.
  - 5 slots are available for: 7.3, 8.5, 10.0, 11.5, 13.0mm
  - Additional stoppers may be pre-installed on the following drill to reduce surgery time and increase convenience.

## Implantation Interval Inspection

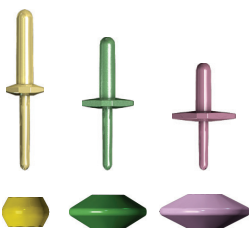

In multiple cases, intervals and directions between implants are measured with MD Parallel Pin. The distance may be measured with in range from  $\varnothing 4.0$  to  $\varnothing 9.0$ , with mesiodistal size of crown determined.

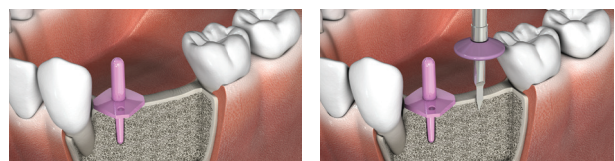

With mesiodistal size of crown determined, Guide Point Drill with MD Spacer installed is used for drilling. The MD Parallel Pin is installed to measure intervals between implants.

# Crestal Approach Tool Components

## • Point Drill



Diameter	Ø1.35
Product name	PD13


## • Initial Drill



Diameter	Ø2.2	Ø2.2
Product name	TSD22CS	TSD22CL

## • MD Parallel Pin

Diameter	Product name
Ø4.0 / Ø5.0	MDPP4050
Ø6.0 / Ø7.0	MDPP6070
Ø8.0 / Ø9.0	MDPP8090

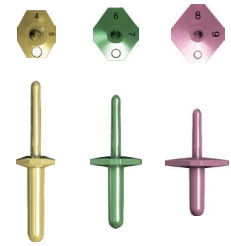

## • MD Spacer




Diameter	Ø4.0	Ø7.0	Ø8.0
Product name	MDS 4025	MDS 7030	MDS 8030

## • Cortical Tap (IS-II active)



Ø3.5	Ø4.0	Ø4.5	Ø5.0	*Optional
CT38S08	CT43S08	CT45S08	CT50S08	

## • Cortical Tap (IS-III / IT-III active)

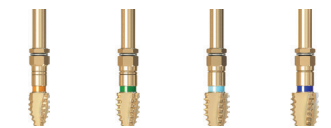

Diameter	Ø3.5	Ø4.0	Ø4.5	Ø5.0
Product name	CT38S09	CT43S09	CT45S09	CT50S09

## • Drill



Diameter	Ø3.0	Ø3.0	Ø3.5	Ø3.5	Ø4.0	Ø4.0	Ø4.5	Ø4.5
Type	Short	Long	Short	Long	Short	Long	Short	Long
Product name	TSD30CS	TSD30CL	TSD35CS	TSD35CL	TSD40CS	TSD40CL	TSD45CS	TSD45CL

## • Direction Pin



Diameter	Ø3.6	Ø4.5
Product name	DP36	DP45

## • Stopper




Stopper	3.0	4.0	5.0	6.0	6.6	7.3	8.5	10.0	11.5	13
Drilling Length	4.2mm	5.2mm	6.2mm	7.2mm	7.8mm	8.5mm	9.7mm	11.2mm	12.7mm	14.2mm
Product name	DS030C	DS040C	DS050C	DS060C	DS066C	DS070C	DS085C	DS100C	DS115C	DS130C

## • IS Fixture Driver



Ratchet (Short)	Ratchet (Long)	Contra angle (Short)	Contra angle (Long)	*Optional
ISFD10R	ISFD15R	ISFD05C	ISFD05CL	

## • IT Fixture Driver

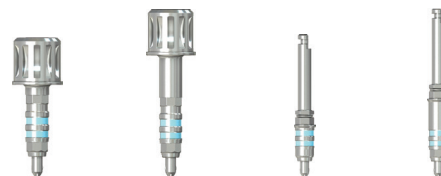

Ratchet (Short)	Ratchet (Long)	Contra angle (Short)	Contra angle (Long)
ITFDO31RS	ITFDO31RL	ITFD31S	ITFDO31L